


VisiClear[®]

Surgical Smoke Evacuator

VisiClear[®] is our most advanced acute care offering for surgical smoke evacuation. Combining Safeport Technology[™], an occlusion warning, filter life tracking, and one-touch programmable procedure modes for open tubing, electrosurgical pencil, and laparoscopic procedures, VisiClear effectively captures and filters surgical smoke, removes odors, particulates and other potentially hazardous by-products.



SYSTEM 5000
GENERATOR

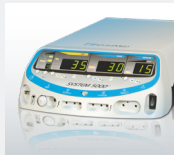

PLUMEPEN[®] ULTRA
PENCILS



EZLINK[®]
AUTOMATION ACTIVATION DEVICE



VIROSAFE[®] FLUID TRAP
FILTER



To learn more about these and other innovative products, call **800-448-6506** or visit **CONMED.com**.

Features & Benefits

Procedure Operating Modes

- One of three operating modes (Laparoscopic, Open Tubing, and Electrosurgical Pencil Mode) can easily be selected, automatically customizing filter life based on surgical procedure.

Advanced Sound Control Design

- Whisper Technology™ ensures ultra-quiet operation with unprecedented sound quality.

Sleek Shelf Design

- Ergonomically designed to accommodate stacking.

Occlusion Warning

- Safety feature that alerts medical personnel when an occlusion has occurred.

Intuitive Color Touch Screen

- Easy to use, requires minimal training, and can be operated with gloved hand.

Unique Sliding Port Closure

- SafePort Technology™ ensures reduced risk of staff coming in contact with a contaminated port.

Technical Information

Noise Levels (Nominal Values)

@ minimum power < 43 dBA

@ maximum power < 55 dBA

Dimensions (H x W x D)

7.1 in x 13.9 in x 16.1 in

18.0 cm x 35.4 cm x 40.9 cm

Flow Rates

30 cfm (839 lpm) using 7/8 in (22 mm) tubing

Weight

With filter installed: 17.5 lbs (7.9 kg)

Without filter installed: 14.5 lbs (6.6 kg)

Static Suction

150 mmHg vacuum

Certifications

UL, cUL, CE

Warranty

Three (3) Year Full Warranty

To learn more about this and other innovative products, call **1-800-448-6506** or visit **CONMED.com**.

Ordering Information

	Description	Qty	Catalog Number
NEW	VisiClear® (100/120 Volt System) Includes one filter	1/Each	VC120
	VisiClear® (200/240 Volt System) Includes one filter	1/Each	VC220
	VisiClear® (100/120 Volt System) Does not include filter	1/Each	VC220 NF
	VisiClear® (200/240 Volt System) Does not include filter	1/Each	VC220 NF

Distributed by:

CONMED Corporation
525 French Road, Utica, NY 13502

Customer Service: 1-800-448-6506
International Sales: 1-315-797-8375
Fax: 1-800-438-3051

CONMED.com
info@mail.CONMED.com