



FERDINAND MEDICAL

Advanced Equipment, Trusted Care

VCare® Plus and VCare® DX Plus


VCare® Plus helps surgeons deliver superior control of the anatomy, effective uterine manipulation and identification of key surgical landmarks for a safe, efficient and repeatable procedure.



AIRSEAL® iFS
INSUFFLATION MANAGEMENT
SYSTEM



AIRSEAL®
ACCESS PORTS
VALVE-FREE ACCESS

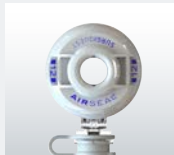

MICROLAP®
INSTRUMENTS
LESS INVASIVE APPROACH



To learn more about these and other innovative products, call 800-448-6506 or visit CONMED.com

VCare® Plus and VCare® DX Plus

New and enhanced VCare® uterine manipulator helps you deliver an even better surgical experience.

- Provides superior exposure and access to pelvic cavity
- Delineates vaginal fornices for easy identification of colpotomy plane
- Maintains pneumoperitoneum throughout the procedure

Guide Ridge

Bright green cup with guide ridge helps define the precise incision landmark for colpotomy

Optimized Distal Cup Design

Enhanced mobility and articulation
Cup geometry and multiple sizes allow for better seating of the cervix
Harmonic compatible

Natural Curved Shaft

Conforms to the patient's anatomy for easier insertion and better exposure

Graduation Marks

Aids in attaining proper depth of insertion

Proximal Occluder

Back cup locks its position to secure the forward cup and prevents loss of pneumoperitoneum

Ergonomic Handle

Improved control for better manipulation and access to key surgical targets

A more comfortable hand position to reduce fatigue

NEW

Easy To Use

One-piece, fully assembled device reduces time required for assembly.
Eliminates cost associated with additional components, repair, sterilization, or lost parts.

Additional Benefits

Allows for suturing of the cup for improved stability of the uterus.
10cc inflatable balloon reduces risk of uterine perforation by anchoring the uterus and stabilizing the shaft within the uterine cavity.




FERDINAND MEDICAL

Advanced Equipment, Trusted Care

VCare® PLUS
AND VCARE® DX PLUS

VCare® DX Plus provides excellent control of the uterus without distention of vaginal fornices. Ideal for:

- Tubal Ligation
- Myomectomies
- Diagnostic Laparoscopy
- Treatment of Endometriosis

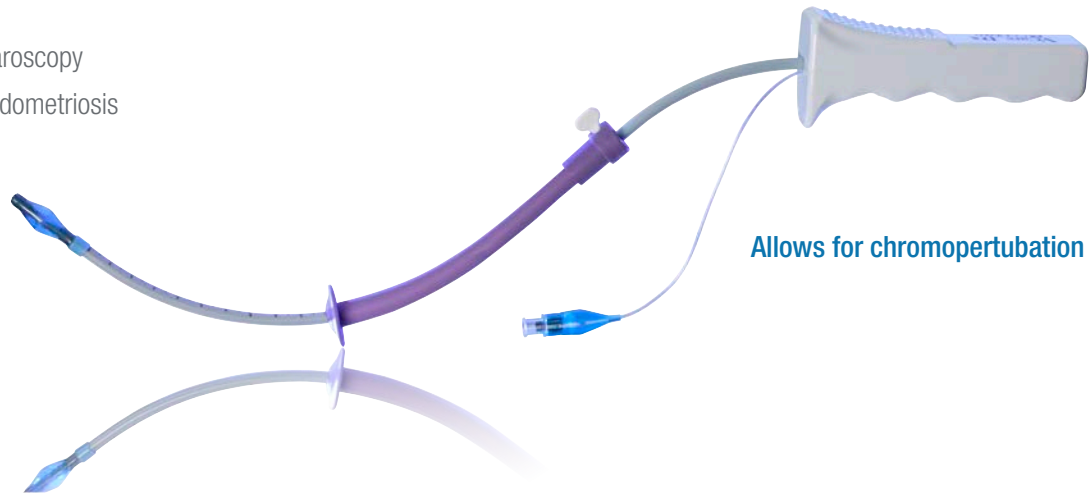

Allows for chromoperturbation

Ordering Information

Description	Qty	Catalog Number
VCare® Plus		
VCare® Plus, Small (32mm) Cervical Cup	8/cs	60-6085-200A
VCare® Plus, Medium (34mm) Cervical Cup	8/cs	60-6085-201A
VCare® Plus, Large (37mm) Cervical Cup	8/cs	60-6085-202A
VCare® Plus, X-Large (40mm) Cervical Cup	8/cs	60-6085-203A
VCare® DX Plus		
VCare® DX Plus Manipulator	8/cs	60-6080-000A

Customer Service: 1-800-448-6506
International Sales: 1-315-797-8375
Fax: 1-800-438-3051

CONMED.com
info@mail.CONMED.com