



FERDINAND MEDICAL

Advanced Equipment, Trusted Care

www.ferdinandmedical.com



Omnia Stapling Platform

Ultraflex 90

90° Articulation Reload



Omnia
Platform



90°
Articulation



Pro Height
Staples



30 / 45 / 60 mm
Staple Length



CE
Certificated



Industry-first
90° Articulation

Witness a monumental pivot in stapling technology.

Offering unparalleled 90-degree articulation, Ultraflex 90 reloads provide surgeons with exceptional freedom of maneuverability, allowing them to reach farther and access difficult areas with unprecedented precision.


Specifications

Product Code	Staple Line Length	Trocar Compatibility (mm)	Color	Staple Height (Open/Closed) (mm)	Tissue Thickness	Number of staple rows	QTY/ BOX
IDSTV30TAN/ IDSTV30TANB	30mm	12	Tan	2.0/2.5/3.0 0.75/1.0/1.25	Vascular/ Thin	6	6
IDSTV45TAN/ IDSTV45TANB	45mm						
IDSTV60TAN/ IDSTV60TANB	60mm						
IDSTV30PUL/ IDSTV30PULB	30mm		Purple	3.0/3.5/4.0 1.25/1.5/1.75	Medium/ Thick		
IDSTV45PUL/ IDSTV45PULB	45mm						
IDSTV60PUL/ IDSTV60PULB	60mm						
IDSTV45BLK	45mm	15	Black	4.0/4.5/5.0 1.75/2.0/2.2	Extra Thick		
IDSTV60BLK	60mm						

Optimized for Lower Anterior Resection (LAR)

The Ultraflex 90 excels in confined, narrow spaces, making it particularly suited for Lower Anterior Resection (LAR) procedures. Surgeons can maintain a safe distance from the tumor site while minimizing the risk of damaging surrounding vital structures.

Lower Reach in Confined Spaces

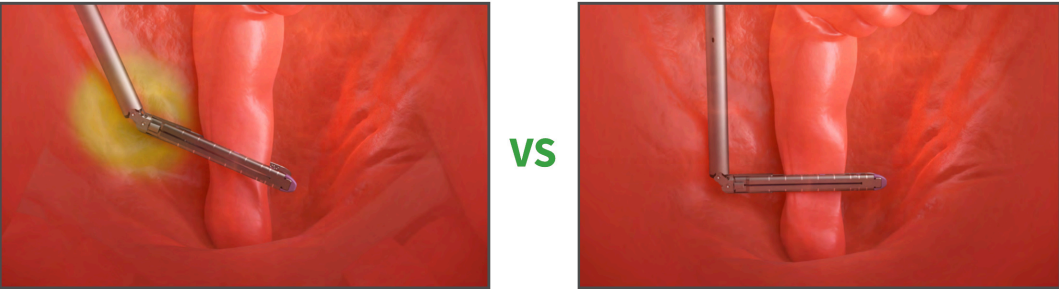


Conventional Reloads

Conventional reloads with 45-degree articulation can be limiting when positioning at lower angles.

Ultraflex 90 Reloads

Ultraflex 90's unique 90-degree articulation allows surgeons to reach deeper into the pelvic region, ensuring a clean transection away from the tumor site with improved accessibility.

Reduced Leakage Risk

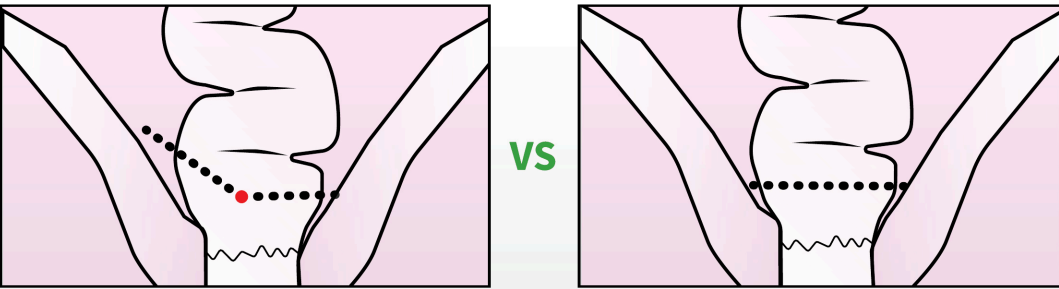

Conventional Reloads

Conventional reloads often require multiple firings, resulting in crossing staple lines that significantly increase the risk of leakage.

Ultraflex 90 Reloads

Ultraflex 90 achieves a clean transection with a single firing, reducing the risk of leakage caused by crossed staple lines.

Flexible applications

The Ultraflex 90 offers surgeons increased flexibility across a wide range of surgeries. Tissue approximation is made easier, with less tension and more access than ever before.



Colorectal



Esophageal



Thoracic

Contacts

Ferdinand Medical ApS
8600 Silkeborg
Denmark

T +45 80 29 10 00
E info@ferdinandmedical.dk
W www.ferdinandmedical.com

Importer



Distributor

