



FERDINAND MEDICAL

Advanced Equipment, Trusted Care

www.ferdinandmedical.com



The Best Choice

Medscope Minimally Invasive Surgery Medical Device



330 mm
Working length



5 and 10 mm
Shaft diameter



10 / 15 / 20
Titanium clips




Remaining-clip
indicator

SOLID CLIP™ CA16 Series | 10mm

Single Use Clip Applier

SOLID CLIP™ applier is a disposable single-use sterile surgical instrument which provides the modern minimally invasive surgery a best occlusion solution


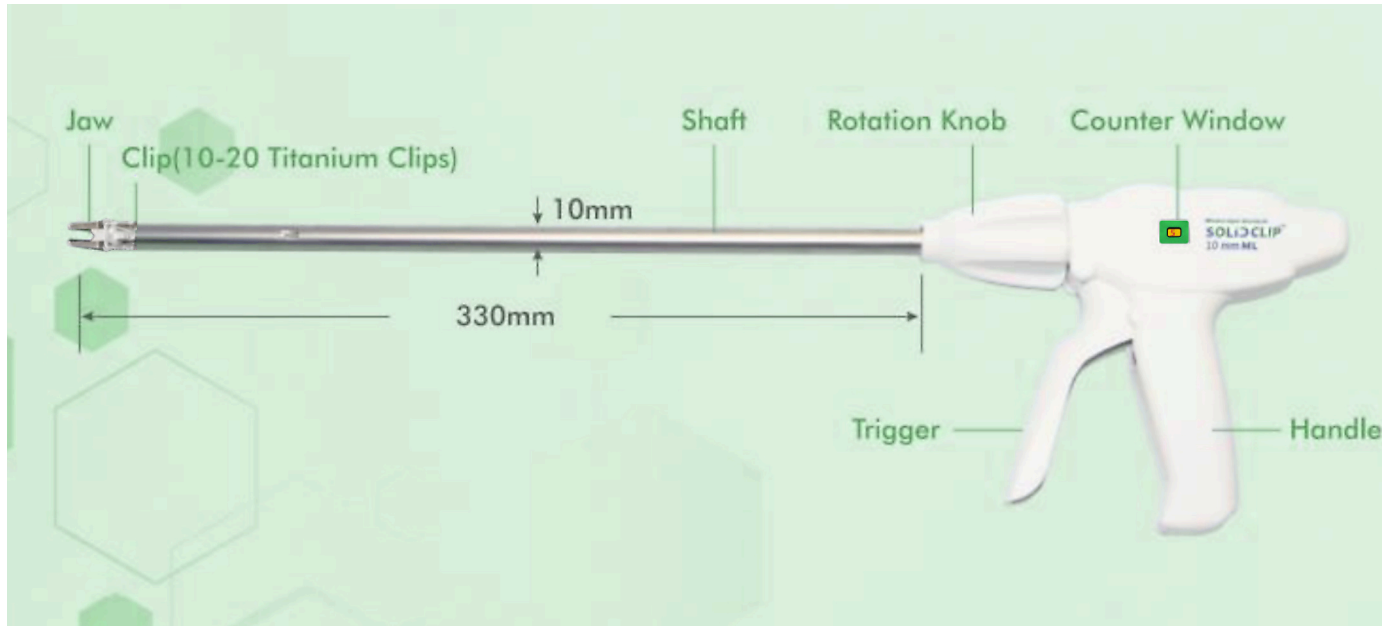


Specification

Model No.	Shaft Diameter	Clip Size	Number of Clips	Packaging
CA165	10mm	Large	10	5pcs/box, 4 boxes/carton
CA166	10mm	Large	15	5pcs/box, 4 boxes/carton
CA167	10mm	Medium/Large	15	5pcs/box, 4 boxes/carton
CA168	10mm	Medium/Large	20	5pcs/box, 4 boxes/carton
CA169	10mm	Medium/Large	10	5pcs/box, 4 boxes/carton


Clip Dimensions

Dimension Unit : mm	Clip Size	
	(ML)	(L)
W	4.5	4.5
H	7.8	9.8
L	9	11

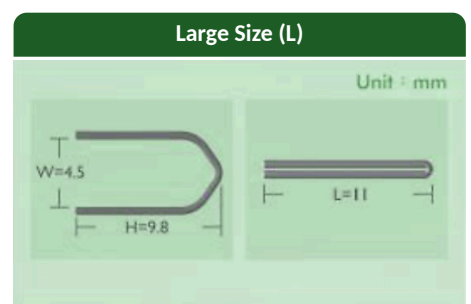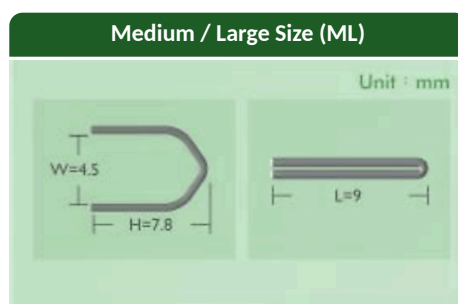

Features

Jaw Design

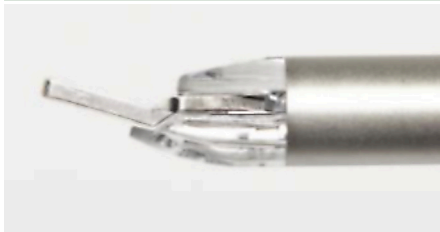

- 15° slop design for better visibility

Jaw Design




- Anti-drop design of the Jaw

Clip Design



- Inner-crossed pattern of the clip offers a firm hold on the vessel

Clip Design

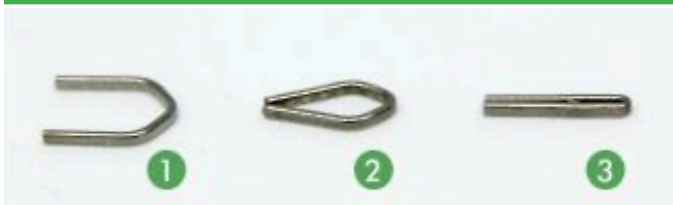

- The closing feature of the clip can provide a proper occlusion of vessel

Clip Material




- Implantable CP Gr.1 titanium clips. The compositions is ASTM F67 Standard Compliant

Clip Counter Window




- Preload 10~20 titanium clips
- The yellow block shown means that there are 5 clips left in this product. "0" means that there are ZERO clips left and the trigger will stop functioning.

Rotation Knob

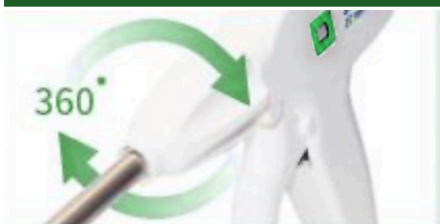

The rotation knob can make the shaft rotating 360° in either direction

Handle

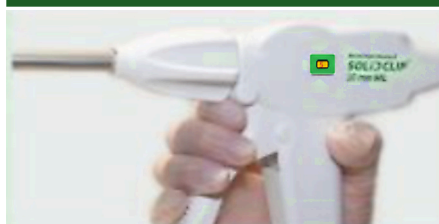

Ergonomic handle design brings excellent control from precise clip loading to closure.

Trigger



- Anti-slip design of the trigger
- Anti-fire design of the trigger when running out of the clips



FERDINAND MEDICAL


Advanced Equipment, Trusted Care

SOLID CLIP™ 5 mm

Single-use Clip Applier



Patented Pitch-Loading design

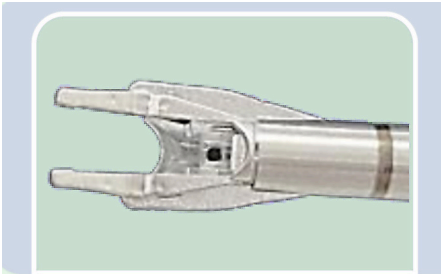

Anti-drop Design

of the
Jaw Structure.




Pitch-Loading
Clip feeding Design

Lighter & Stable
Ligating Force.



Interlock Design

when running out
of the clips.

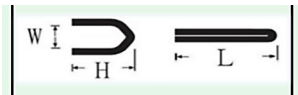

Model	Clip Size	W	H (mm)	L	Preloaded Clip Number	Color Indicator			Dimensions
						Clip Number			
									
CA-568	M/L	3.5	8.2	9.0	15 Clips	1st ~ 12th	13th ~ 15th	No Clips	5 x 330 mm

SOLiDCLIP™



Medscope Laparoscopy

Contacts

Ferdinand Medical ApS
8600 Silkeborg
Denmark

+45 80 29 10 00
info@ferdinandmedical.dk
www.ferdinandmedical.com

Importer



Distributor

